

TRITEC™

Contact Layer Dressings Featuring AFM® Technology



Every Milliken TRITEC and TRITEC Silver dressing utilizes Active Fluid Management® Technology created by Milliken Healthcare Products, LLC. This three-part design combines two layers of high-performance fabric using a micro-knit process, creating a moisture transfer mechanism that pulls away and traps excess exudate. The combination of these fabrics keeps the wound bed moist, the periwound skin healthy, and harmful exudate away from a healing wound.


TRITEC Silver combines the moisture management benefits of AFM Technology's three-part design with a proprietary ceramic silver ion technology. This powerful, controlled release of ionic silver inhibits a broad spectrum of infection-causing microbes including antibiotic-resistant strains such as MRSA. The release of silver into the dressing is activated by the body's own sodium and begins inactivating pathogens within 24 hours of dressing application. Antimicrobial activity is indicated for up to seven days.

INDICATIONS FOR USE

- Partial thickness burns
- Incisions
- Skin grafts
- Donor sites
- Lacerations
- Abrasions
- Stage I-IV dermal ulcers (vascular, venous, pressure and diabetic)

BENEFITS

- Removes excess exudate through Active Fluid Management Technology
- Minimizes risk of maceration
- Protects healthy, intact skin
- Allows for fewer dressing changes
- Improves patient comfort through conformability

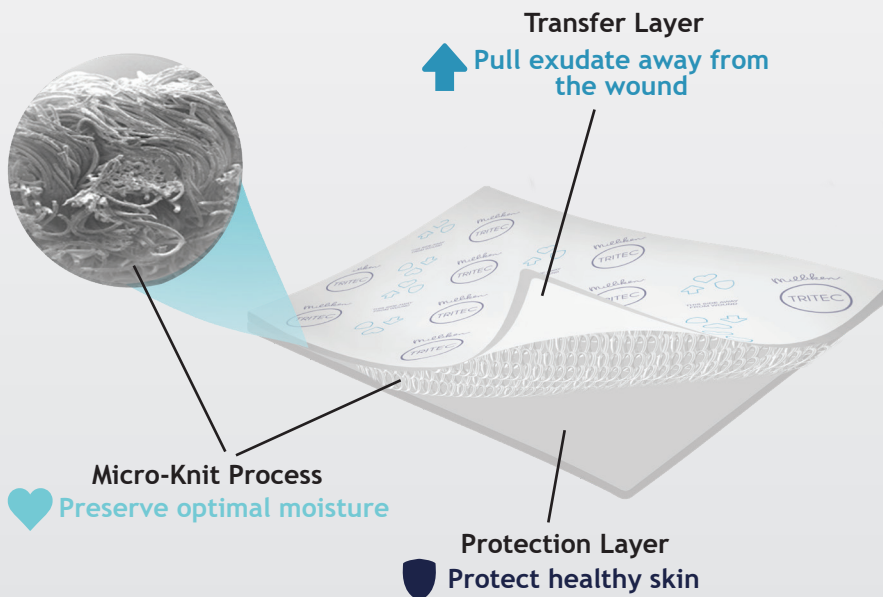

Milliken™

MILLIKEN HEALTHCARE PRODUCTS, LLC

ORDERING INFORMATION




TRITEC	Size	Code	HCPCS	Packaging
	4" x 5"	3000020040	A6207	10 ea/bx 10 bx/cs
	4" x 5" Pre-Cut	3000112816	TBD	10 ea/bx 10 bx/cs
	4" x 48"	3000020051	A6208	10 ea/cs



TRITEC Silver	Size	Code	HCPCS	Packaging
	1" x 24"	3000004574	A6207	5 ea/bx 10 bx/cs
	4" x 5"	3000005137	A6207	10 ea/bx 10 bx/cs
	4" x 5" Pre-Cut	3000112331	TBD	10 ea/bx 10 bx/cs
	6" x 6"	3000004573	A6207	10 ea/bx 10 bx/cs
	4" x 48"	3000004575	A6208	10 ea/cs
	8" x 16"	3000005211	TBD	50 ea/cs

Three Easy Ways to Order:

1. Preferred Distributor
2. Preferred DME
3. Direct from Milliken

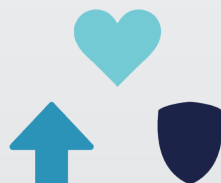

pull preserve protect

Pull exudate away from the wound.
Preserve optimal moisture.
Protect healthy skin.

medproducts@milliken.com | healthcare.milliken.com | 1-866-491-6556