ASSIST Silver Absorbent

Antimicrobial Low-Adherent Absorbent Dressing Featuring AFM® Technology




Milliken,

MILLIKEN HEALTHCARE PRODUCTS, LLC

Every Milliken ASSIST Silver Absorbent dressing utilizes Active Fluid Management Technology[®] created by Milliken Healthcare Products, LLC. This three-part design combines two layers of high-performance fabric using a micro-knit process, creating a moisture transfer mechanism that pulls away and traps excess exudate in the two highly absorbent layers incorporated into the dressing. The grey impervious film backing prohibits moisture from leaking through the back of the dressing, but can easily be removed to allow for airflow and / or use of additional secondary dressings.

ASSIST Silver Absorbent combines the moisture management benefits of AFM's three-part design with a proprietary ceramic silver ion technology. This controlled release of ionic silver inhibits a broad spectrum of infection-causing microbes including antibiotic-resistant strains such as MRSA. The release of silver into the dressing is activated by the body's own sodium with antimicrobial activity for up to seven days.

The patented low adherent layer in ASSIST Silver Absorbent eliminates the need for a separate non-adherent layer, which saves O.R. and nursing time during application and maintenance.

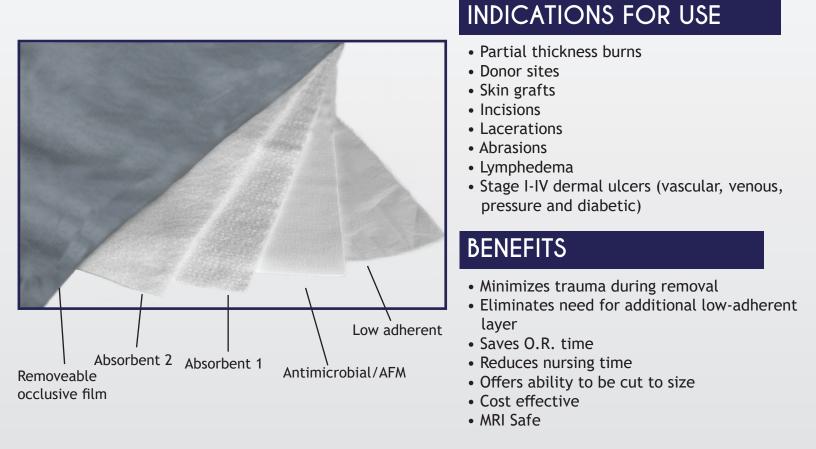



ORDERING INFORMATION

Size	Code	Packaging
26" x 26" Extremity	3000117008	15 ea/cs
36" x 24" Vest	3000117007	10 ea/cs
36" x 48" Half Bed	3000117009	5 ea/cs
18" x 24"	3000116362	40 ea/cs
24" x 36"	3000006964	20 ea/ca

Three Easy Ways to Order:

- 1. Preferred Distributor
 - 2. Preferred DME
- 3. Direct from Milliken

medproducts@milliken.com | healthcare.milliken.com | 1-866-491-6556